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Dear customer,

Thank you for choosing our **pore cleanser**. If you have any questions, please contact customer service via our website, at:

www.dspro.de/kundenservice

Explanation of the symbols

Danger symbol: This symbol indicates danger of physical injury. Please read and follow the safety instructions carefully.

Symbol for direct current

Explanation of the signal words

WARNING warns of *possible* serious injuries and danger to life

CAUTION warns of mild to moderately severe injuries

NOTICE warns of material damage

Intended use

- The product is intended for pore cleansing (removal of sebum and blackheads), microdermabrasion and to stimulate blood circulation in the skin. The skin is sucked in by the product. Due to the vacuum, sebum and impurities are extracted from the pores and blood circulation in the skin circulation is stimulated.
- Warning, do not misuse! Do not use the product to clean ears.
- The product is designed for private use, not for commercial use and not for medical or therapeutic purposes.
- Use the product only for the specified purpose and only as described in these instructions. Any other use is regarded as improper use.

SAFETY INSTRUCTIONS

- **WARNING – danger of suffocation!** Keep the packaging away from children and animals.
- **WARNING – danger of chemical burn from battery acid!** Batteries can be life-threatening if swallowed. Keep batteries and product out of reach of children and animals. If a battery has been swallowed, seek medical help immediately.
- **CAUTION – danger of chemical burn from battery acid!** If a battery leaks avoid physical contact with the battery acid. Wear protective gloves if necessary. Should this occur, rinse the affected areas immediately with plenty of clean water, and consult a doctor immediately.
- **CAUTION – risk of injury!**
 - Keep the product away from children.
 - Do not use the product for rosacea, birthmarks, warts, injured, diseased, (sun) burned or very sensitive skin.
 - Do not use in case of illness / dental disease and wounds.
 - People with insensitive skin should not use the product.
 - The following people should use the product only after medical consultation: People with a pacemaker or an electronic implant, others who are in a risk group or have concerns about their health.
 - When using, avoid the eye area and the thyroid area.
 - **Move the product constantly during use!** Holding it in place for a few seconds may result in a vacuum bruise (“love bite”).
 - If an unpleasant sensation occurs while using the product, stop using it immediately!
 - If the skin colour on the hand holding the product changes, stop using it immediately! Prolonged use could lead to hand-arm vibration syndrome.

NOTICE – risk of material and property damage

- There must be no liquid in the product. Do not immerse in water! Never use when bathing or showering!
- Remove a leaking battery from the product immediately. Do not use the product if a battery is leaking.
- Only use the battery type specified in the technical data.
- Remove batteries from product if they are used or if the product will not be used for a long time. This helps prevent damage from leaks.

- Always replace all batteries. Do not use different battery types, brands or batteries with different capacities.
- Only use original accessories from the manufacturer so as not to impair the functioning of the product and to prevent possible damage.
- Protect the product from heat, open flames, direct sunlight, sub-zero temperatures, prolonged humidity, moisture and impacts.
- Do not leave the product switched on unattended.
- Keep the product away from flammable materials and substances.
- In case of damage / malfunction, switch off the product immediately. Do not use the product if it or an attachment is damaged.
- Do not disassemble the product or try to repair it yourself. In the case of damage, contact customer service.


Included in delivery

1 x main unit	4 x silicone rings, big
1 x round pore attachment, large	2 x silicone rings, small
1 x microcrystalline attachment	5 x filters
1 x oval attachment	1 x operating instructions
1 x round pore attachment, small	

2 batteries (1.5 V, type AA) are also needed. These are not included in the delivery.

Check the scope of delivery for completeness and transport damage. In the event of damage do not use; contact customer service

Product overview

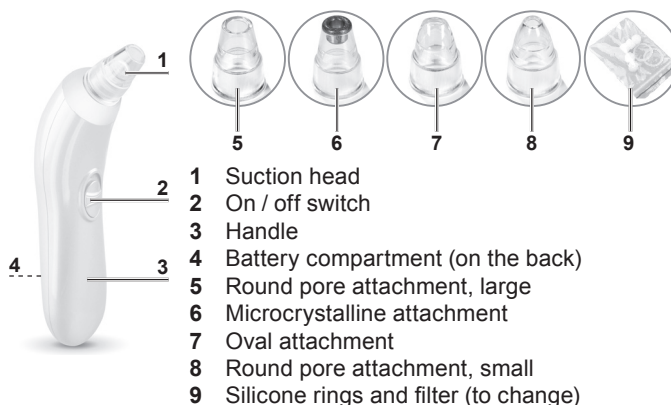Insert / change batteries

Before inserting the battery, check whether the contacts in the product and on the battery are clean, and clean them if needed.

1. To open the battery compartment: Slide open the battery compartment cover in the direction of the arrow and remove it.
2. Remove used batteries if necessary and replace them with two new batteries (1.5V; type AA). Make sure the poles are aligned correctly (+/-)!
3. Put on the battery compartment cover and slide it shut in the opposite direction to the arrow.

Replace the batteries if power decreases.

Usage

PLEASE NOTE

- For hygiene reasons, the product should only be used by one person. Thoroughly clean the attachment used immediately after use.
- Do not use peeling products on affected areas of the skin before or after use.
- Do not use the product more than once a week on dry skin and twice a week on oily or combination skin.
- Do not use for more than 5 minutes.
- Test the product on the forearm before first use. If the skin does not react undesirably, the product can be used on the face.
- Tingling and reddening of the skin is normal during use and should disappear within half an hour.
- **To prevent blemishes forming, keep the attachment steady on the skin during use.** It is recommended to stretch the skin, e.g. on the cheeks, with your fingers, so it is not sucked in.
- Move the product slowly in strips. **Do not circle and do not press the skin!**

Remove make-up before treatment and clean the face as usual. Clean the face after each treatment and apply a care product.

Four attachments are included in delivery:

- **Large round pore attachment** for pore cleansing, not suitable for the eye area
- **Small round pore attachment** for pore cleansing in places that are difficult to reach, e.g. in the nose area, not suitable for the eye area
- **Oval attachment** to stimulate blood circulation in the skin
- **Microcrystalline attachment** for microdermabrasion: Removal of dead skin flakes by friction and suction on the face, elbows, knees, etc.

Pore cleansing

Tip: For best results, pores should be opened before use. Dip a small clean towel in hot water, wring it out and put it on the face. Alternatively, you can carefully hold the face over a bowl of hot water or use the product directly after showering or bathing.

1. Dab the face dry.
2. Put the round pore attachment on the suction head.
3. Slide the on / off switch upwards. The product switches on.
4. To remove **sebum and blackheads**, place the pore attachment on the skin and move it slowly over the skin in small strips from the middle of the face to the outside. Only pull the strips in one direction. If necessary, repeat the process. Briefly suck in **acne**, then quickly remove the pore attachment again. Do not move to the sides!
5. After treatment, switch off the product by sliding the on / off switch downwards.
6. Remove the pore attachment and clean it.

Microdermabrasion

1. Dry the face well.
2. Put the microcrystalline attachment on the suction head.
3. Slide the on / off switch upwards. The product switches on.
4. Start with facial treatment on the chin and work your way to the forehead. Slowly move the microcrystalline attachment outward in strips from the middle of the face. Only pull the strips in one direction and do not press the product! **Please note:** The microcrystalline attachment can also be used on body parts such as elbows or knees.
5. After treatment, switch off the product by sliding the on / off switch downwards.
6. Remove the microcrystalline attachment and clean.

Stimulate skin circulation

1. Dry the face well.
2. Put the oval attachment on the suction head.
3. Slide the on / off switch upwards. The product switches on.
4. Briefly suck in skin areas where wrinkles form (e.g. around the mouth), then quickly remove the oval attachment. Do not repeat too often in one place, so that blemishes don't form!
5. After treatment, switch off the product by sliding the on / off switch downwards.
6. Remove the oval attachment and clean.

Cleaning and storage

PLEASE NOTE

- Never use corrosive or abrasive cleaners.
- Clean the attachments under running water immediately after each use. If necessary, use a mild detergent and cotton buds.
- If necessary, wipe the product with a damp cloth and then dry it.
- Keep the product away from dust, dirt and direct sunlight in a dry place away from children and animals.

Maintenance

- If the attachments on the suction head do not hold, one of the silicone rings may be damaged and must be replaced. Delivery includes replacement silicone rings.
- Replace the filter in the suction head if it is dirty. Delivery includes replacement filters.

Change filter

Similar to the one shown.

1. Pull the insert out of the suction head (see image).
2. Use a pointed object (e.g. needle, tweezers) to pry the dirty filter out of the suction head.
3. Insert a fresh filter into the suction head.
4. Insert the insert back into the suction head.

Technical data

Article number:	09215
Model number:	GWF174
Power supply:	3 V DC (2 x 1.5V; type AA)
Operating instructions ID:	Z 09215 M DS V1 0419

Disposal

Dispose of packaging materials in an environmentally sound way and take recyclables to collection facilities.

Remove batteries before disposing of the product and take them to the appropriate collection site separately from the product.

This product is subject to European Directive 2012/19/EU on Waste Electrical and Electronic Equipment (WEEE). Do not dispose of the product as normal household waste, but through an officially approved disposal business according to environmental regulations.

Customer Service / Importer:

DS Produkte GmbH, Am Heisterbusch 1, 19258 Gallin, Germany • ☎ +49 38851 314650 *)

*) Calls to German landlines are subject to your provider's charges.

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All rights reserved.
DS Produkte GmbH, Am Heisterbusch 1, 19258 Gallin, Duitsland • ☎ +49 38851 314650 *)